

SciLampPowerControl User Manual

Table of Contents

1. Introduction	3
2. Installation.....	3
a. Software Installation Sequence	3
b. Uninstalling Sciencetech Lamp Power Control:	6
3. Software Description	7

List of Figures

Figure 1. Lamp power control screen with labelled features. Please not that not all the functionality comes wit the default software. Depending on the accessories purchased, more functionalities are enabled.	7
Figure 2. Setup Power Supply screen	8

1. Introduction


The Lamp Power Control is an application which allow the user to control a Sciencetech Power Supply.

2. Installation


Sciencetech Lamp Power Control is Windows 7, 10, and 11 compatible. It may work on Windows 8 but is not thoroughly tested in that operating system.

a. Software Installation Sequence


1. You can run the installer exe from a network location, USB or from the computer's hard drive. If you copy the installer to a hard drive before installation makes sure you copy the folder structure exactly. For example, in the image on the right you must copy all 3 items to the same directory before you can run
Sciencetech_LampPowerControl.exe

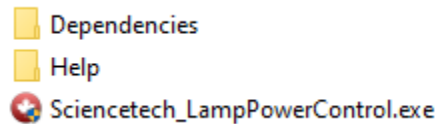

2. Launch
"Sciencetech_LampPowerControl.exe".
In order to run, the installer will request administrator privilege. Based on your system, you will see a popup window asking for approval, please click on "Run"(Windows 7) or "Yes" (Windows 10).

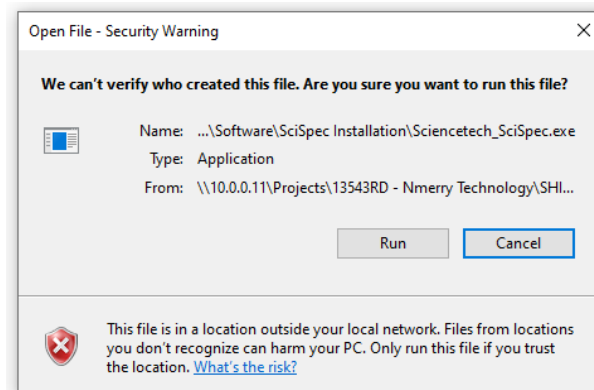


3. Review the license agreement.
Click 'I Agree' to continue

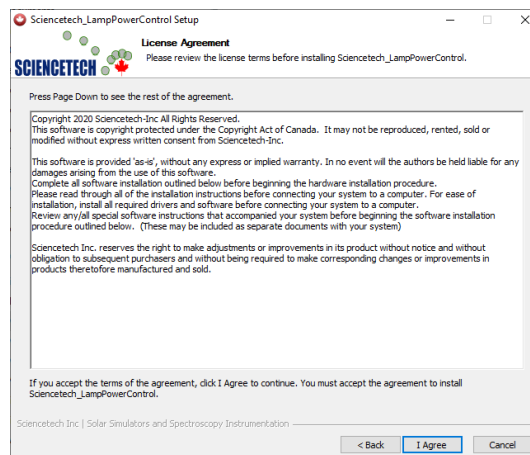


4. Choose to create a StartMenu Shortcut. This will create a shortcut for SciMono application and SciMono User manual.

Note: At any time, a shortcut will be created on your desktop.

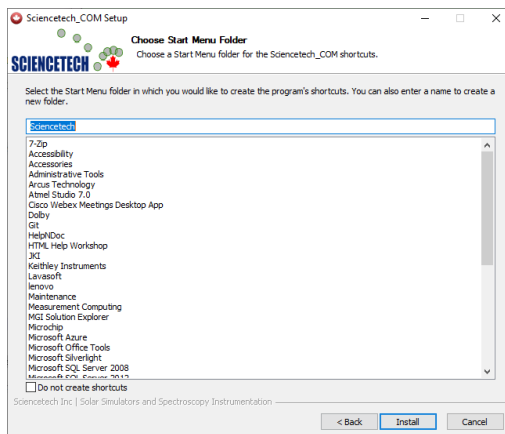


5. Continue to install the software:
If existing configuration files are detected, the installer will give you two options:

- **To not overwrite** any existing file, in that case, it is possible the old configuration file may not be compatible with this installation.
- **Overwrite the existing file** (Recommended)
The previous one will be saved in a backup folder on your desktop.

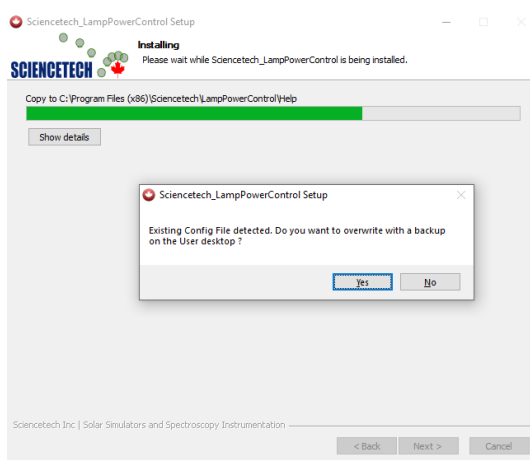


6. Install third party drivers
Click on extract and then follow the instruction (click on next and accept the agreement)

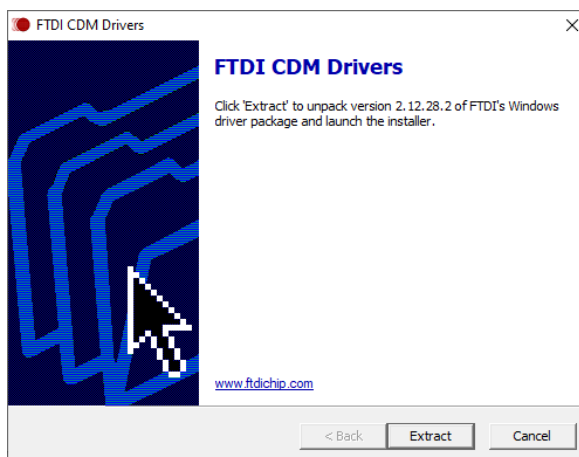


7. Successful Software Installation

When the software completed a successful installation, you will see the final popup

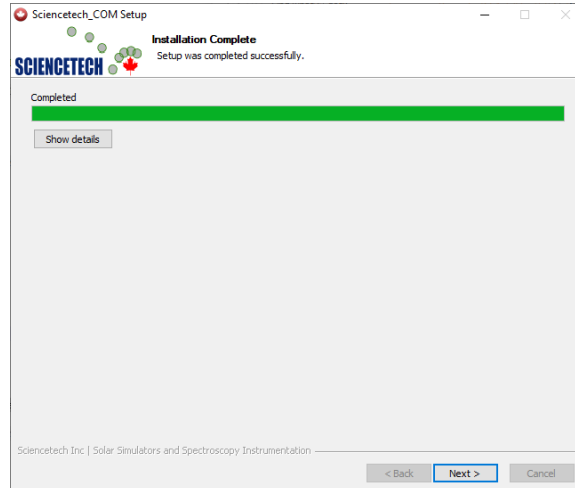

8. At closure, you can start your new software

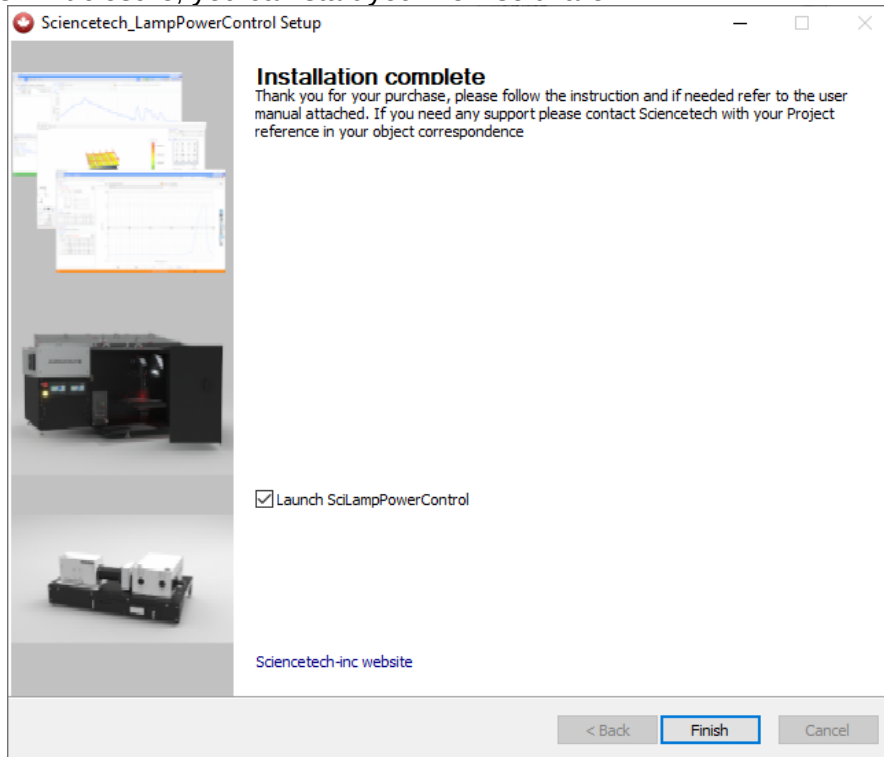


9. User Manuals:


By default user manuals will be copied to the following directory: (C:\Program Files (x86)\Sciencetech\LampPowerControl\Help)

b. Uninstalling Sciencetech Lamp Power Control:


1) If you want to uninstall Lamp Power control, please go to the Add/Remove program.
Look for Lamp Power Control in the list of installed programs.

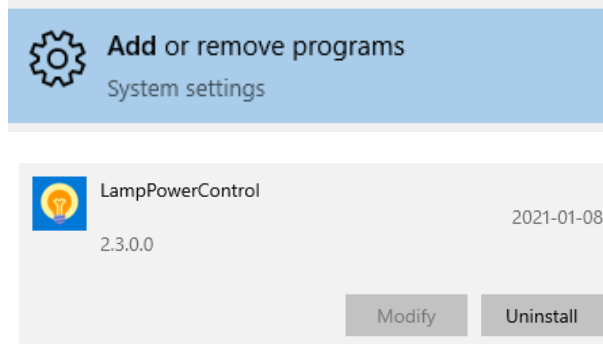

2) Click on uninstall and follow the instructions.

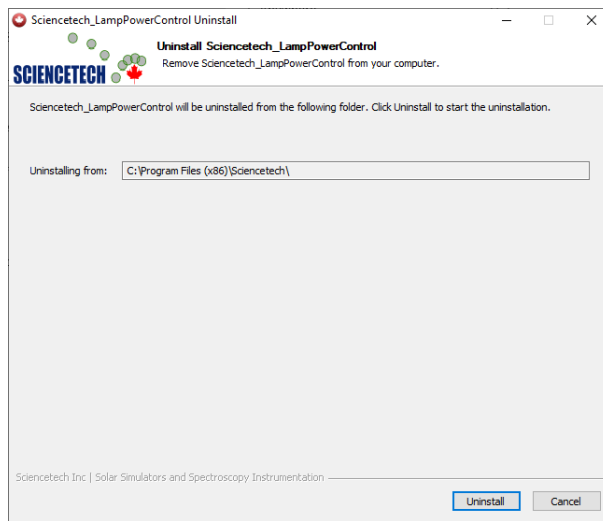

Note, if you have requested the uninstallation of the config file, a popup will ask you for confirmation and if you want to keep a backup on the desktop. (Recommended if you have modified the initial config file provided with the installation).

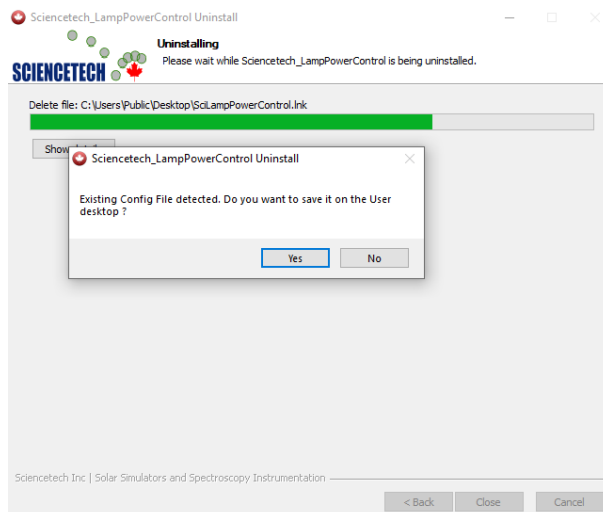

3. Software Description

1. All connections from the power supply to the lamp should be verified. Electrical hookup is detailed in your solar simulator or lamp house user manual. Connect the touchscreen power supply RS232 to your computer with the provided cable. Determine the COM port number for the touchscreen power supply. COM port information can be accessed through Device Manager on your Windows computer.
2. Double click the SciLampPower.exe shortcut to run.
3. Open the settings window by clicking the settings button and enter the correct COM port for the power supply connected to your computer in the LampPowerCommPort text box and click the Save button.
4. To initialize the power supply, click the Initialize button.
5. The operation of the GUI follows the same operation as the touchscreen on the power supply. Please see the accompanying Lamp House or Solar Simulator user manual. See the following diagram for a list of the button functions.

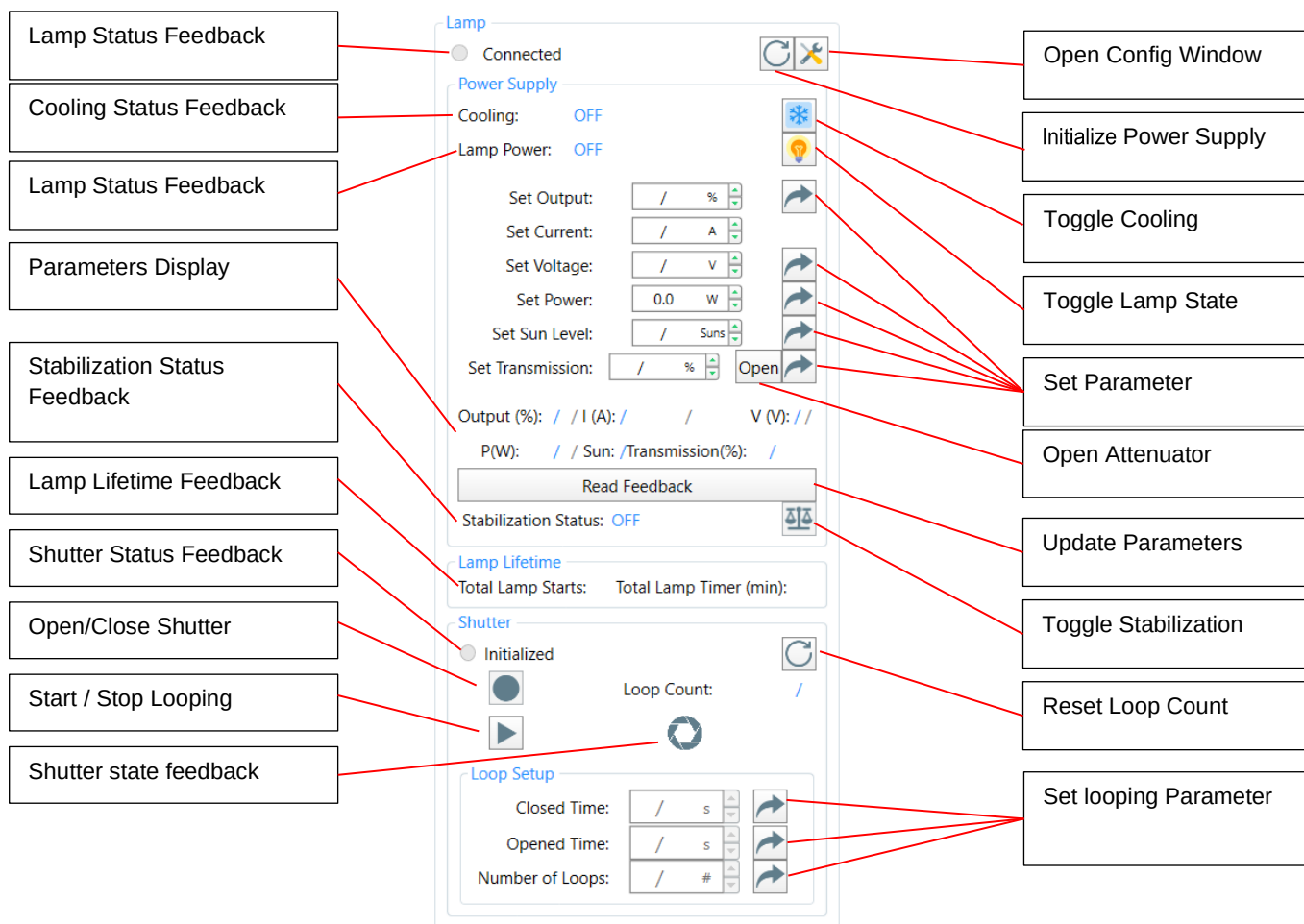


Figure 1. Lamp power control screen with labelled features. Please not that not all the functionality comes wit the default software. Depending on the accessories purchased, more functionalities are enabled.

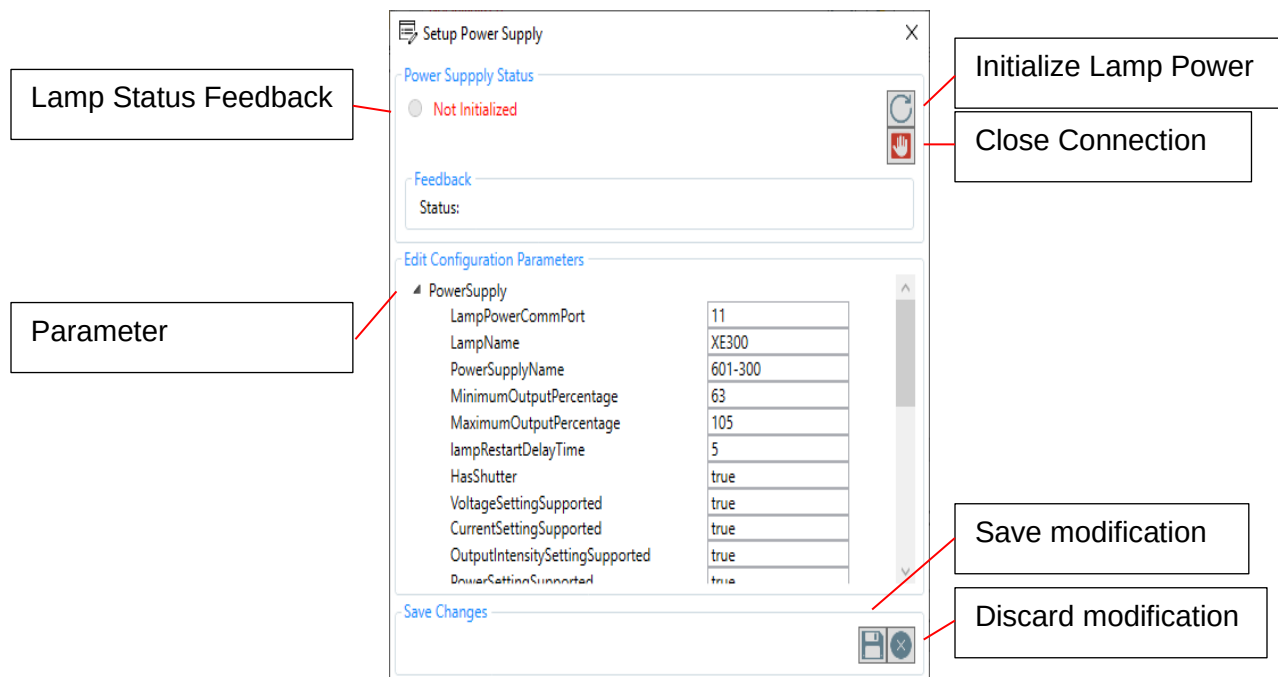


Figure 2. Setup Power Supply screen

WARNING: Do not change MaximumVoltage without first contacting Sciencetech.